

OLYMPUS[®]



EVIS EXERA DUODENOVideoscope

OLYMPUS TjF TYPE 160F

SETTING THE STANDARD IN DUODENOSCOPY



With an extra-large 4.2 mm diameter channel
to accommodate wider stents
and improved angulation handling,
this versatile duodenoscope
is ideal for a variety of applications




Powerful Features Offer Proven Results

- Superb treatment capability with exceptionally large 4.2 mm diameter channel that accommodates a wide range of Endo-Therapy accessories including a variety of stents.
- New distal-end configuration puts the objective lens and forceps elevator closer together for improved operability in cannulation.
- Newly designed angulation knobs easy to turn with light handling power for better operability.
- Slim 11.3 mm insertion tube.
- Ergonomically designed grip enhances scope maneuverability while redesigned easy-to-access knobs and user-selectable switches improve operability.
- Fully compatible with the CV-160/140/100.
- Scope ID function stores individual scope information and displays it on the monitor to facilitate endoscopy suite management.





Specifications

Optical System	Field of view	100°
	Direction of view	5° Backward oblique-viewing
	Depth of field	5 to 60 mm
Distal End	Outer diameter	13.2 mm
	Outer diameter	11.3 mm
Bending Section	Angulation range	Up 120°, Down 90°, Right 110°, Left 90°
Working Length		1240 mm
Total Length		1550 mm
Instrument Channel	Inner diameter	4.2 mm
	Minimum visible distance	10 mm
Endo-Therapy accessory entrance/exit position in field of view		

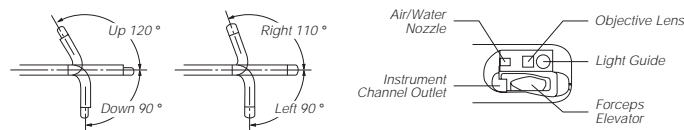

Olympus business areas

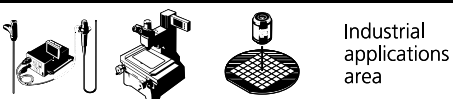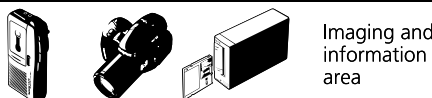

OLYMPUS
FOCUS ON LIFE

<http://www.olympus.com>

OLYMPUS OPTICAL CO., LTD.

San-Ei Building, 22-2, Nishi Shinjuku 1-chome, Shinjuku-ku, Tokyo, Japan

OLYMPUS OPTICAL CO. (EUROPA) GMBH
Postfach 10 49 08, 20034 Hamburg / Wendensstrasse 14-16, 20097 Hamburg, Germany

OLYMPUS AMERICA INC.
2 Corporate Center Drive, Melville, N.Y. 11747-3157, U.S.A.

OLYMPUS KEYMED
Keymed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England

OLYMPUS SINGAPORE PTE LTD.
491B River Valley Road #12-01/04, Valley Point Office Tower, Singapore 248373

OLYMPUS BEIJING REPRESENTATIVE OFFICE
Rm. 818, South Tower, Beijing Kerry Centre, No.1 Guanghua Road, Chaoyang District, Beijing, 100020, China

OLYMPUS MOSCOW LIMITED LIABILITY COMPANY
117071, Moscow, Malaya Kaluzhskaya 19, bld. 1, fl.2, Russia

OLYMPUS AUSTRALIA PTY. LTD.
1/104 Ferntree Gully Road, Oakleigh, VIC 3166, Australia